


AIR2CARE

LINET

A Great Deal of Pressure Relief

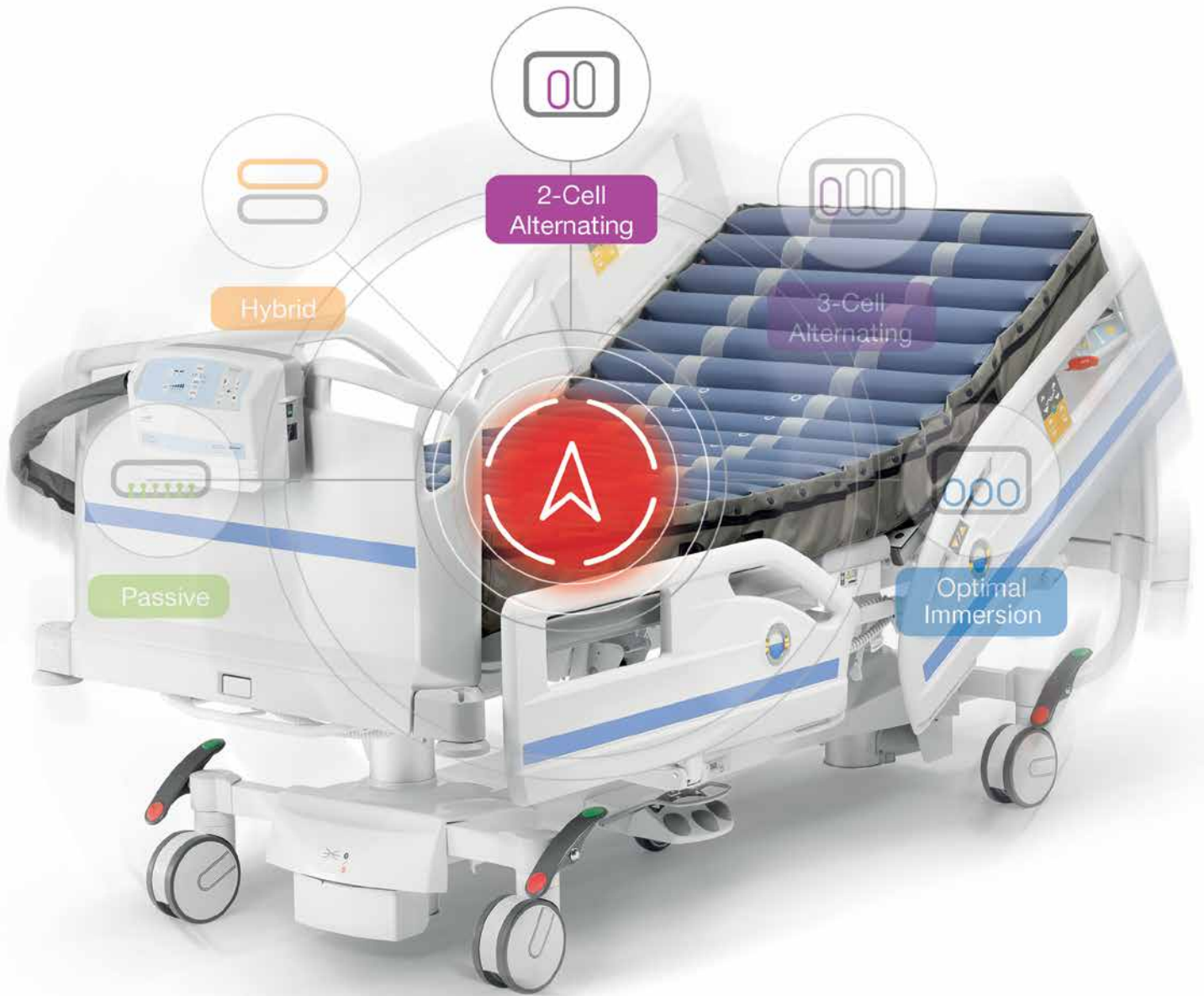


2-cell alternating pressure mattress family

AIR2CARE

A Great Deal of Pressure Relief

EFFICIENCY IN THE PREVENTION OF PRESSURE INJURIES

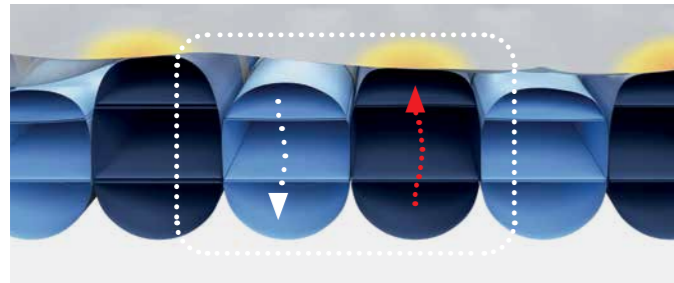

2-cell alternating system provides intermittent pressure reduction.



Fowler Boost helps to minimize the risk of bottoming out.

SIMPLICITY OF HANDLING BROUGHT BY USER-FRIENDLY FEATURES



Easy Smart System with automatic pressure setting adjusting to the patients weight and position.




Easy patient transportation supported by Transport mode that can be enabled by simple disconnection from System Control Unit(SCU).



Convenient maintenance provided by individual cells that are removable and can be easily cleaned, decontaminated or replaced.

SAFETY AND COMFORT AT THE CORE OF PATIENT CARE

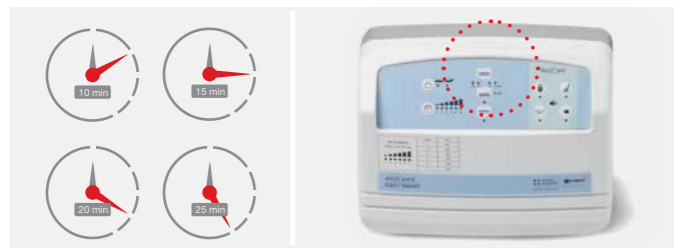

Wide choice of cycle times based on patient's condition and individual preference.



Integrated compressor provides less noise and vibration for better patient comfort.

INTEGRATED MATTRESS

The compressor of the mattress is integrated into the Essenza bedframe, along with an easy to use control panel available on the bed siderails.



Integrated solution



MATTRESS REPLACEMENT

Alternating Air2Care mattress replacement with an external compressor can be placed in any bed frames with compatible dimensions.



Mattress replacement



TECHNICAL SPECIFICATION

MODELS AIR2CARE 8 AIR2CARE 6 AIR2CARE 5 AIR2CARE 4

FEATURES

Air Decks	Air+Air	Air+Air	Overlay	Overlay Seating Cushion
Modes	Alternate, Max, Constant Low Pressure	Alternate, Max, Constant Low Pressure	Alternate, Max, Constant Low Pressure	Alternate, Max, Constant Low Pressure
Transport mode	Yes	Yes	Yes	Yes
Cycle Time*	10, 15, 20, 25 min	10, 15, 20, 25 min	10, 15, 20, 25 min	10, 15, 20 min
Fowler Boost	Single point	Single point	Single point	—

MATTRESS DIMENSION

Mattress replacement (cm)	200 x 86 x 20	200 x 86 x 15	200 x 86 x 12,5	48 x 46 x 10
	200 x 100 x 20	200 x 90 x 15		
	200 x 120 x 20	210 x 86 x 15		
	200 x 140 x 20	210 x 90 x 15		
	210 x 86 x 20	220 x 86 x 15		
Integrated mattress (cm)	200 x 86 x 20	200 x 86 x 15		

ADDITIONAL PARAMETERS

Maximum patient weight (kg)	86 cm wide: 180 kg	160	140	130
	100 cm wide: 250 kg*			
	120 cm wide: 280 kg*			
	140 cm wide: 300 kg*			
Minimum patient weight (kg)	40	40	40	40
Cover	Waterproof and Vapor Permeable, Stretch Cover	Waterproof and Vapor Permeable, Stretch Cover	Waterproof and Vapor Permeable, Stretch Cover	Waterproof and Vapor Permeable, Stretch Cover
Zipper	360°	360°	360°	270°
Handles	Yes	Yes	Yes	—

* Related to mattress replacement only



Members of LINET Group

LINET spol. s r.o.

Želevčice 5 | 274 01 Slaný | Czech Republic

tel.: +420 312 576 400 | fax: +420 312 522 668 | e-mail: info@linet.com | www.linet.com



www.linet.com