

HYBRIMATT

LINET

The Hybrid Path to Efficiency



Active and passive systems in one

WHY HYBRID MATTRESS?



Hybrid mattress systems are support surfaces that contain elements of both, a passive and an active system.

Hybrid mattresses allow caregivers to easily step up or step down as the patients risk of developing pressure injuries changes throughout their hospital journey. It's versatility ensures that patients can remain on a single mattress.

HybriMatt contributes to pressure injury prevention, streamlines hospital operations, and offers a variety of mattress options, making it


the Hybrid Path to Efficiency.



HYBRID APPROACH FOR PRESSURE INJURY PREVENTION



Active system with 2-cell alternation provides periodic pressure reduction.

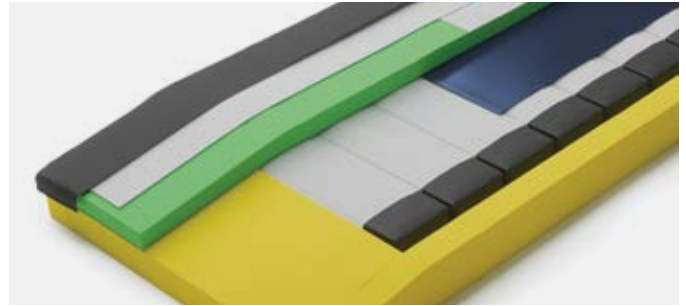

Passive system helps to redistribute pressure thanks to a viscoelastic foam layer in the mattress.

ENHANCING EFFICIENCY IN HOSPITAL PROCEDURES



Hybrid mattresses **simplify daily procedures for caregivers** and support flexibility in prevention strategies.



Easy **mattress cleaning and inspection** provided by full access to all the components of the mattress.

VERSATILE MATTRESS OPTIONS



Integrated mattress option provides a cohesive solution for Essenza bed.



Mattress replacement can be placed on any bedframe.

TECHNICAL SPECIFICATION

MODELS



HYBRIMATT 100



HYBRIMATT 200



HYBRIMATT 100 IN



FEATURES

| | | |
|---------------------------------|--|---|
| Air Decks | Air+Foam | Air+Foam |
| Foam Layer | 5 cm viscoelastic foam 10 cm polyurethane foam | 5 cm viscoelastic foam 10 cm polyurethane foam |
| Modes | Alternating Pressuere Therapy, Max, Constant Low Pressure | Alternating Pressuere Therapy, Constant Low Pressure |
| Cycle time mattress replacement | 10, 15, 20, 25 min | 12 min |
| Cycle time integrated mattress | 10 min | |
| Microclimate Maintenance | No | Yes |

MATTRESS DIMENSION

| | | |
|---------------------------|--|--|
| Mattress replacement (cm) | 200 × 86 × 15; 200 × 90 × 15; 208 × 86 × 15 | 200 × 86 × 15; 200 × 90 × 15; 208 × 86 × 15 |
| Integrated mattress (cm) | 200 × 86 × 15 | |

ADDITIONAL PARAMETERS

| | | |
|--|---|--|
| Maximum patient weight (kg) | 230 | 230 |
| Minimum patient weight (kg) ¹ | 40 | 40 |
| Cover | Comfort cover, Endurance cover ² | Comfort cover, Endurance cover, Zoned cover |
| Zipper | 270°, 360° ² | 360° |
| Handles | Yes | Yes |
| Evacuation cover (option) | – | Yes |
| X-ray sleeve (option) | – | Yes |

¹ Related to mattress replacement only. ² Related to integrated solution only.

