

Solido II

Universal bed server



Patient

- A more pleasant bed stay
- Tilt option for reading



Staff

- Easy clean concept
- Simple handling



Management


- Modern design
- Ergonomic movement by gas spring



Solido II

The Solido II server is a practical accessory allowing patients greater comfort when eating and spending their leisure time, such as writing or reading. The benefit of the server is its universal use with a bed or chair.

- 1 DOUBLE-SIDED**
tiltable eating board
(optional)
- 2 HEIGHT SETTING**
using the control handle
- 3 DOUBLE CASTORS**
Tente 50 mm or 75 mm
(optional)
- 4 BRAKES**
on 2 castors

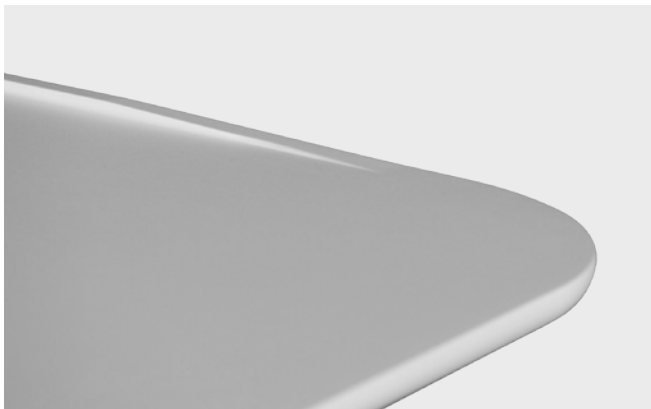

DESIGN

The patient and nurse will appreciate the modern look and simple control. The eating board has clean lines, The server frame is either silver RAL 9006 or white RAL 9002.

SIMPLE CONTROL

The Solido II server can be equipped with a spring for upward direction control without the use of the control handle.




PRACTICAL TILT

The server is equipped on both sides with a tiltable eating board which patients will appreciate when reading or watching a tablet. The tiltable board can be locked into five positions: -40 , -20 , 0 , $+20$ and $+40$ degrees. The server can also be supplied in a version with a fixed eating board without tilting.



PROTECTIVE CATCH

prevents the sliding and falling of a book or other objects.

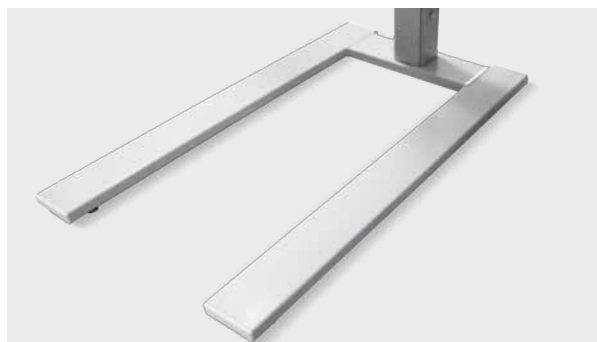

EXTRA LOW CLEARANCE

The Solido II server option with integrated castors offers extra low clearance.

Technical data

Solido II	Non-tiltable	Tiltable
Eating board height adjustment (+5/-10 mm)		
50 mm castors	683-983 mm	713-1013 mm
75 mm castors	700-1000 mm	730-1030 mm
Undercarriage height		
50 mm castors	102.4 mm (62.4 mm between floor and undercarriage)	
75 mm castors	119.4 mm (79.4 mm between floor and undercarriage)	
Undercarriage width	420 mm	420 mm
Eating board dimensions	450 x 960 mm	450 x 992 mm
Maximum weight on eating board	30 kg	30 kg
Tilting angle	N/A	-40°, -20°, 0°, 20°, 40°
Handle position Solido II	left sided (optional right sided)	

Low clearance version	Non-tiltable	Tiltable
Height adjustment of eating board (+5/-10 mm)	700-1000 mm	730-1030 mm
Undercarriage height	34.3 mm (14.3 mm between floor and undercarriage)	
Undercarriage dimensions	390 x 810 mm	390 x 810 mm
Eating board dimensions	450 x 960 mm	450 x 992 mm
Maximal weight on eating board	30 kg	30 kg
Tilting angle	N/A	-40°, -20°, 0°, 20°, 40°
Handle position	universal	



Members of LINET Group

LINET spol. s r. o.

Želečnice 5 | 274 01 Slaný | Czech Republic
 tel.: +420 312 576 400 | fax: +420 312 522 668 | e-mail: info@linet.com | www.linet.com

